

HemoCue® Hb 201+ System



The Standard for Uncompromising Accuracy

Providing lab accuracy and ease of use, the HemoCue® Hb 201+ system has become a standard in Hb point-of-care testing. Healthcare providers around the world rely on the immediate results so they can make the right decisions when they need them most – right at the point of care.

With dedicated support and service as well as unmatched training and education based on 30 years of experience, you can count on HemoCue for the right solutions for all your needs.

Accuracy Starts With Us



To learn more about the HemoCue® Hb 201+ system, please scan the QR code with your smartphone or visit hemocue.com

Have Confidence in Your Answers at the Point of Care

- ▶ Precise factory calibration against the ICSH reference method
- ▶ Patented microcuvette technology with excellent lot-to-lot reproducibility
- ▶ Built-in selftest
- ▶ Blood-based liquid controls available


Get Easy Access to Lab-Quality Accuracy

- ▶ Capillary, venous or arterial whole blood sample
- ▶ Brief training with virtually no maintenance
- ▶ Link result with patient ID for medical record integration
- ▶ Printer interface




HemoCue® Hb 201+ System

Principle	Modified azidemethemoglobin reaction; dual wavelengths (570 nm and 880 nm) for compensation of turbidity
Calibration	Factory calibrated against the ICSH reference method; needs no further calibration
Sample Material	Capillary, venous or arterial whole blood
Measurement Range	0–25.6 g/dL (0–256 g/L, 0–15.9 mmol/L)
Results	Within 60 seconds
Sample Volume	~10 µL
Dimensions	160×43×85 mm (6.30×1.69×3.35 inches)
Weight	350 g (0.77 pounds) with batteries installed
Storage Temp.	Analyzer: 0–50 °C (32–122 °F) Microcuvettes: unopened 15–30 °C (59–86 °F); three-month open vial stability
Operating Temp.	15–30 °C (59–86 °F)
Power	AC Adapter or 4 AA batteries
Interface	Printer and HemoCue® Basic Connect™ including optional barcode scanner
Quality Control	Built-in selftest ; optional liquid controls
CPT Code	85018
Product Item Number	Hb 201+ Analyzer 121721 Hb 201 Microcuvettes (200/box) 111716 Hb 201 Individually Packaged Microcuvettes (100/box) 111715

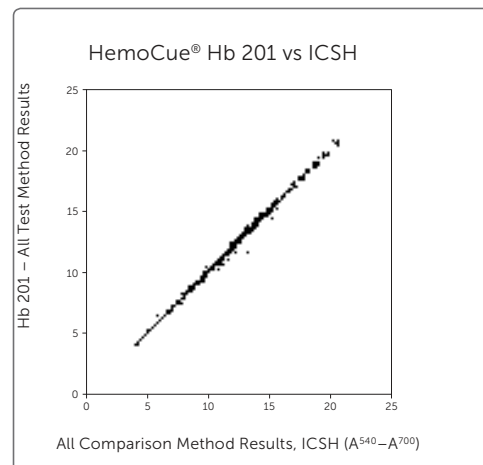


The CPT codes provided are based on AMA guidelines and are for informational purposes only. CPT coding is the sole responsibility of the billing party. Please direct any questions regarding coding to the payer being billed.


Three Simple Steps

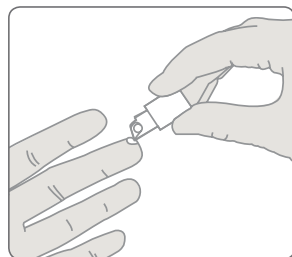

1 Fill microcuvette.

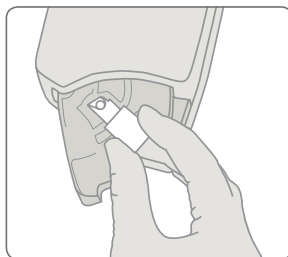

2 Place microcuvette into analyzer.

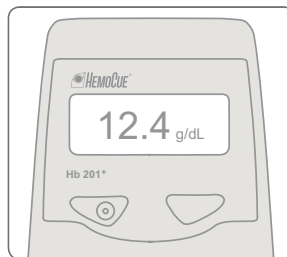

3 View results.

HemoCue America | California, USA
Phone 800.881.1611 | Fax 800.333.7043 | info@hemocue.com | hemocue.com

HemoCue has been a leader in point-of-care medical diagnostics for over 30 years. We specialize in giving healthcare providers lab-quality accuracy with results comparable to that of a clinical lab.

